

TRIANGO 100 Duo

■ ■ ■ ■ ■ ■ ■

medicinelaser[®]
LATAM



reddot winner 2020

The reliable double luminaire system. TRIANGO 100 DUO

Two treatment fields? Low room heights? With the TRIANGO 100 duo ceiling luminaire, the premium treatment light from Derungs, various challenges can be solved.

Two light heads, each with 100,000 lux per meter and excellent lighting technology make it the perfect lighting solution for demanding users.



TRIANGO 100 Duo	technical data
luminous power *	2 × 100,000 Lux / 1.0 m
colour temperature *	fixed: 4300 K adjustable: 3700 K / 4300 K / 4700 K
lighting field @ 1.0 m *	Ø fixed 18 cm focusable: 18 – 28 cm
depth of illumination *	80 cm
colour rendering index	Ra > 95
usage	touch panel on cardanic system
class of protection	I
protection rating	IP 43 (luminaire head)



special features	
TRIANGO 100-1 C Duo	dimnable
TRIANGO Fokus 100-3 C DUO	• dimmable • 3 colour temperatures • focusable
fastening	ceiling

TRIANGO 100 Duo



First class light

The specially developed diamond optics relieve the eyes and ensures the best colour homogeneity. Each luminaire head provides 100,000 lux at one meter. Touch pad controls allow the illumination intensity can be adjusted. Optionally the illumination field can be focused or different colour temperatures (3700 K / 4300 K / 4700 K) can be used.

High mobility

The high-quality support arm system offers perfect mobility: an exceptional radius of action of > 180 cm, joints with 360° rotation and a large angle of movement. The aluminum arm system is especially durable and easy to clean.

Product highlights




Duo ceiling light
enables interventions on 2 treatment fields.
Can also be used at low room heights.



Unique Duo arm system



Parallel bearing arm connection
perfect mobility with large radius of action.

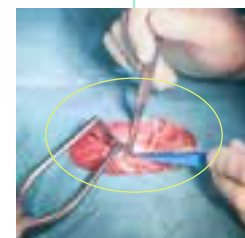

Homogeneous and low-shadow
Illumination field
fixed 18 cm, optional
focusable 18 - 28 cm.

medicinelaser
LATAM
www.medicinelaser.com



+ OFICINAS CENTRALES

ARGENTINA
Av. Doctor Honorio
Pueyrredón 707, Caballito,
CABA, Argentina.
M: +54 11 3404 7158

CHILE
Sucre 1355, Ñuñoa,
Santiago, Chile.
O: +56 2 2269 1981
M: +56 9 4409 7105

ECUADOR
Av. República 1058,
Edificio La Fontana,
Quito, Ecuador.
M: +593 99 120 4597

PERÚ
Calle Almirante Lord
Cochrane 322,
Miraflores, Lima, Perú.
O: + 51 1 704 80 14
M: +51 983 552 813

ventas@medicinelaser.com